

Smile Line

made in switzerland



TRANSITIONS

...4 years passed incredibly fast since we have published the last edition of Smile Line's general catalogue ! 4 years during which we could realize that our society is now facing huge challenges and exceptional changes on many aspects.

It is for us an immense happiness to get back to you today with a new catalogue that gathers all the novelties developed and launched during the past few years.

Since almost 30 years, Smile Line has been considered by the industry a pioneer and a leader for highly exclusive products, for their quality and innovative aspect. 2023 will quickly show to be a milestone thanks to the new and revolutionary concept Optishade (digital shade measurement of teeth colors - featured on p. 44 to 47) that required years of development.

Yes, digital color measurement of teeth is today one of the great challenges in our branch and we are incredibly proud and happy for having succeeded, with our Partners from Styleitaliano, in setting new references at such high level with Optishade.

We are very excited to accompany you in this discovery via the workshops that we organize in our education center. (Information regularly updated on www.smileline.ch)

wishing you pleasant reading and with all our best from Switzerland!

Olivier Schaeren
CEO & Founder



Education center
on www.smileline.ch

« With Optishade,
we stand far beyond a simple
transition; we are talking here
of an unbelievable revolution
that opens a window onto
a completely new and
fascinating world! »



I CREATE...
WE CREATE



INDEX

INDEX

Optishade Styleitaliano	4-7
Smile Lite MDP2	8-11
Smile Lite	12
K-Lite	13
Flexipalette & Photo CAD	14-15
Sidekick, Implant drivers for dentists	16-17
VeneerME	18-19
Holding Gel	19
My Shade Guide for composites	20-21
Sampler	21
Compo-Vibes	22
Compo-Brush	23
Organ.ez	24
N-Architect, Anterior silicone moulds	25
Golden Section	26-27
Styleitaliano	28

OPTISHADE

DIGITAL SHADE MATCHING:
DIGITAL COLOUR IN THE PALM
OF YOUR HANDS !



SMILELINE 2023

For as long as anyone can remember, developers and passionate dental professionals have been trying to find alternatives to the use of commercial dental shade guides. The best solution so far has been to create individualized samples and shade guides out of the actual materials used, either ceramic or composite. In the present, most of dental shade matching is still analogue and needs to be backed up with a lot of experience from the practitioner for a reliable interpretation.

OPTISHADE STYLEITALIANO provides a new solution for shade matching that surpasses all the other digital color measurement devices on the market with its accuracy and repeatability !

The measured shade is communicated in L*a*b* coordinates, universal values in the world of colorimetry, very intuitive and easily interpretable by the user.

These L*a*b* coordinates are either :


- interpreted by the user
- used as a comparison with the commercial shade guides databases in the OPTISHADE STYLEITALIANO.
- or exported to the Matisse software, which will calculate for you all the recipes and ceramic powders mixings necessary for the given case

The L*a*b* in the dental field SIMPLIFIED FOR BEGINNERS:

Almost ANY tooth color that you can see is simply a mix of grey **Luminosity (L*)**, **red (a*)** and **yellow(b*)**. Optishade measures these 3 values [coordinates] with high precision and repeatability.

READING COLOR COORDINATES



The L*a*b* color space describes numerically all colors in three dimensions. L* is for Luminosity a* and b* for the color-axis green-red and blue-yellow, respectively.

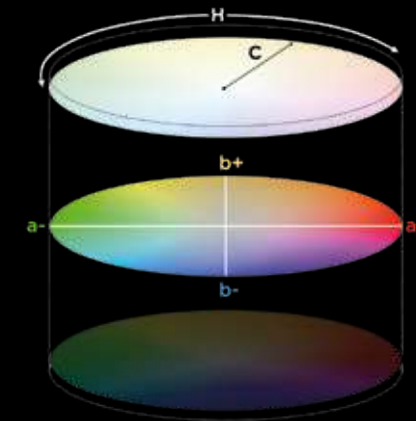

L* scale: Light vs. dark, 100 is white and 0 is black.


b* scale: Yellow vs. blue, where a positive number indicates yellow, and a negative number indicates blue.

a* scale: Red vs. green, where a positive number indicates red, and a negative number indicates green.

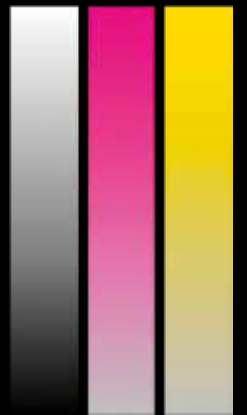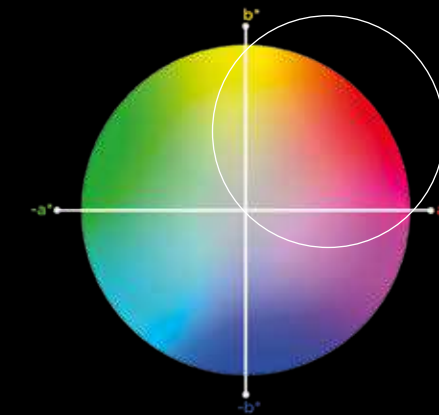

In the dental field, color coordinates are extremely easy to read, as we virtually only have to focus on Luminosity (L*), red (a*) and yellow (b*). It is very uncommon to "visit" the negative numbers.

The L*a*b* scale rises above language barriers enabling anyone to easily communicate color independently of the device, and it provides the necessary information to calculate color differences.

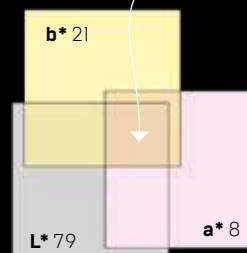


THE DIFFERENT STEPS OF A REAL CASE




Optishade color measurement of the substrate



Optishade color measurement of the adjacent natural tooth

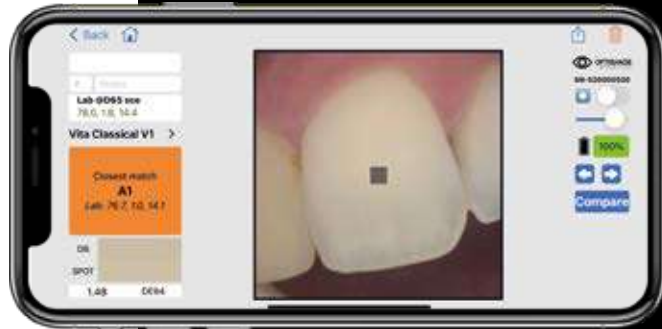

Compare function on the Optishade app : on the right the natural tooth, on the left the crown

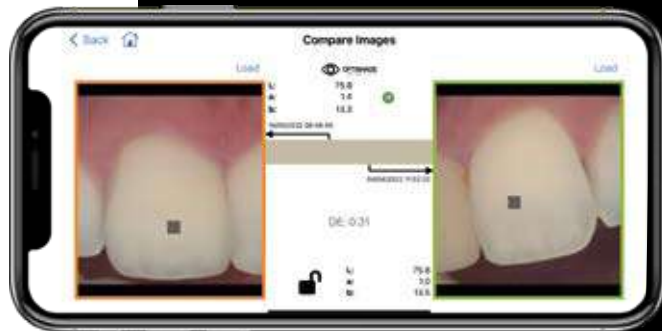

Final result in one attempt. The dental technician never saw the patient



Many thanks to Dr. Marie Clément and Lionel Marslen MDT for the courtesy!

OPTISHADE SET



Description	N° Art.
Optishade Full set	70100-OS
Optishade Capture Guide Cone (5pcs)	70110-OS
Optishade Calibration Cap (5pcs)	70120-OS
Optishade Hardcase (optional)	70130-OS
Optishade Usb-c cable (for ipad pro) (optional)	70140-OS

Only compatible with:



FREE download of your Optishade App (Appstore)

OPTISHADE + MATISSE THE WINNING TEAM !



MATISSE is a software embedded with comprehensive data and artificial intelligence technology to make dental shade matching easy.

You can import the color files from OPTISHADE STYLEITALIANO and Matisse will translate them into simple **ceramic recipes**; will provide you step-by-step solutions for any aesthetic restoration: Aesthetic model-Framework advice- dentin recipes-internal effects-skin enamel.

MATISSE also includes a **FULL MONOLITHIC STAINING STUDIO**, to achieve a perfect shade match extremely easy!

This feature creates staining recipes that you can use for micro-layering, full monolithic crowns and to make color corrections after the first bake with the internal staining technique.



Discover Matisse

SMILE LITE MDP2

MOBILE DENTAL PHOTOGRAPHY

New era for dental photography!

SMILELINE 2023

Smile Lite MDP2 is a plug & play device that is very easy to use and that doesn't need any special App, calibration process or specific education. It will allow you to drastically improve your results with the camera of your smartphone... very easily.

Soon it will become a mandatory tool for you :

- photos for shade taking
- every day communication with the lab
- short video clips
- patients' documentation
- artistic pictures




Smile Lite MDP2 can be considered a "mini photo studio". It is equipped with three groups of latest generation « sunlike » LEDs (light emitting diodes) : each group can be lit up individually, while you also have the possibility to set the power of illumination thanks to a dimmer (four different steps of power).

The back of Smile Lite MDP is equipped with a universal handle / adapter that is adjustable for any model of smartphone having a width between 55-85mm. For matching all tastes and preferences, you will also receive a second phone adapter with no handle.



Universal adapter / handle
55-85 mm width

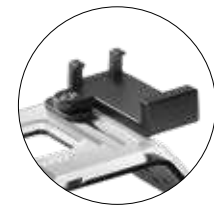

SMILE LITE MDP2 ACCESSORIES




Light Diffuser
Great for softening light and making surface texture more prominent.




Polarizing Filter
Removes light reflection from the teeth and helps to determine the shade and value. It also helps to see the characterizations inside the tooth.

The different lightings


MDP2 with central LEDs at full power, lateral LEDs off

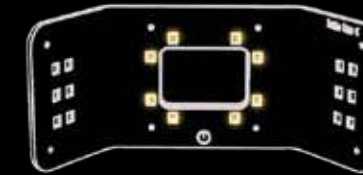

MDP2 with lateral LEDs at full power, central LEDs off

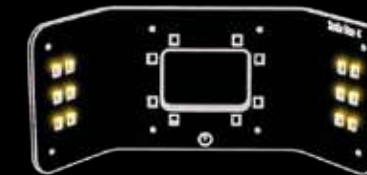

MDP2 with lateral LEDs at full power + diffuser, central LEDs off

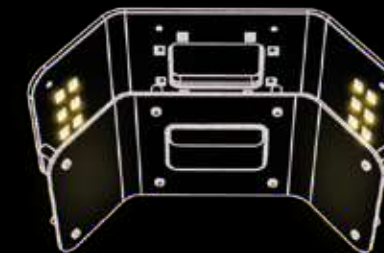

MDP2 with central LEDs at full power + polarizing filter, lateral LEDs off



...and the results



MESSAGE OF THE INVENTOR

Dear colleagues, dear friends,

After some 5 years of experiments with the use of smartphones in dental photography, I had the opportunity to start the development of Smile Lite MDP2 with the community Styleitaliano and the Swiss company Smile Line SA.

It is a huge pleasure and pride for me to introduce you a tool that development was 100% focused on dental photography with smartphones (MDP2 = mobile dental photography). Smile Lite MDP2 is a total revolution for dental photography because of its extreme simplicity of use, the professional results it is possible to reach and also because of its economical aspect. When making the comparison with photography with an SLR (single lens reflex) camera, the Smile Lite MDP2 can be described as a **mini photo-studio** that is able to cover the major needs of the user for:

- A ring flash
- A dual flash
- Photography with softbox
- Photography with a polarizing filter
- Photography using softbox and polarizing filter combined
- Short videos

On the previous page, you can discover a few examples of images that highlight the amazing possibilities when using a smartphone with Smile Lite MDP2. After a very short time of tests and training, you too will be able to achieve such highly professional pictures.

Today on my opinion dental photographs are mandatory in every dental office, for the communication with a patient and with the lab and of course for your own patients' documentation! Smile Lite MDP2 has reached such a level of simplicity that it would be a pity to do without it.

I invite you to be curious, to enjoy playing with light and shadow and discover the phenomenal possibilities of the cameras on today's smartphones.

Wishing you a lot of pleasure at your first steps with Smile Lite MDP2!

Sincerely,

Prof. Louis Hardan DDS, CES, DEA, PhD
 Department of Restorative and Esthetic Dentistry
 Saint-Joseph University, Beirut, Lebanon



SMILE LITE MDP2 SETS



Description	Item no.
Smile Lite MDP2 Master Set	6700-Master



Description	Item no.
Smile Lite MDP2 Basic Set Without Polarizing Filter & Light Diffuser	6700-Basic



MDP2 micro website

SMILE LITE

SHADE MATCHING FOR THE SPECIALIST

Without the proper tool, shade-taking can become very risky and cause a source of failures when producing an aesthetic restoration. Several factors can influence an incorrect perception of true colors, like ambient light (intensity, quality) and the working environment. Even human factors like the level of stress or fatigue on the eye can cause shade-taking errors.

Smile Lite is a revolutionary tool which brings you reliability, simplicity and efficiency – it allows you to **drastically reduce the risk of mistakes during shade-taking.**

Equipped with 5500°K LEDs (light emitting diodes), Smile Lite provides natural and neutral light, stable and reliable regardless of the time of the day or the weather outside. A Li-Ion battery is integrated meaning that the device is easily rechargeable thanks to the mini-USB cable (cable and charger supplied in the set).

Light quality of Smile Lite allows you to **determine with ease and security the hue, value and chroma of natural teeth.**

The rectangular view finder through which the operator observes the teeth clearly limits the “working area” in order to avoid any outside interference.

This “magical” filter amazingly annihilates light reflection (specular and diffuse) and allows the user to **observe the teeth in a totally new way:** easier appreciation of color, better understanding of the depth and transparencies, enhancement of the tiniest details and characterizations.

Observe teeth as you have never seen before!

SMILELINE 2023



Without polarizing filter

With Style Lense polarizing filter



Description	Item no.
Smile Lite, complete set incl. Style Lense polarizing filter	6500
Polarizing filter Style Lense	6510

K-LITE

WHY K-LITE ?

Modern dental composites and adhesives contain fluorophores which help in producing better aesthetic results. This gives them a fluorescent property.


Under certain wavelengths of illumination (UV-A 395 to 400nm), these materials can be distinguished from tooth tissue. Therefore, significantly reducing risk of iatrogenic tooth damage when they need to be removed. The K-Lite is a ground breaking dental device that allows:

1. the removal of excess or unwanted dental materials such as composites or adhesives. Some of the examples include orthodontic brackets, attachments and retainers, residual cements etc.
2. its use in dental trauma (Splint) and/or in the so-called Fluorescence-assisted Identification Technique (FIT). These two objectives were the motivation for the author to develop and produce the new lamp K-Lite in cooperation with the company Smile Line SA, Switzerland.
3. It can assist in the identification of certain bacterial by products found in plaque, calculus and in infected dentine.

- Easy and fast use
- Dual UV and Daylight LEDs
- Portable - No cable
- Easy to handle - Light weight


Make the invisible visible

Description	Item no.
K-Lite, Complete set	3000-KL
K-Lite light guide	3010-KL
K-Lite disposable protective sleeves, 200 pcs	3011-KL

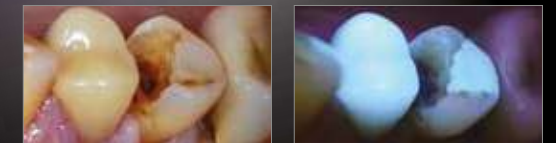


Without UV light it's not easy to detect the kind of restoration material and its extension.

With K-Lite's UV light the identification of composite materials and their extensions is easy.

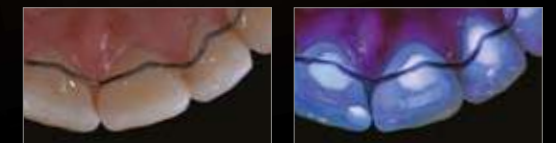

Without K-Lite: no possible composite adhesive detection.

With K-Lite: visibility of the material for control of excess.

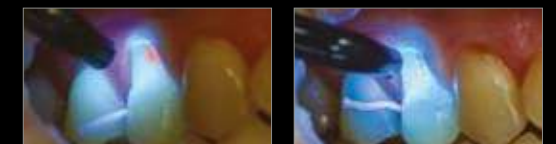

Presence of pigments due to bacteria degradation in a cervical lesion.

Control after removal of infected area.

3000-KL
K-Lite



FLEXIPALETTE

In partnership with Styleitaliano, a study group of passionate clinicians and researchers, Smile Line introduces on the market an assortment of very unique and innovative contrastors intended for intra-oral photography: Flexipalette.

Flexipalette is made out of a bendable copper sheet 100% overmoulded with a medical grade silicone rubber, latex free, safe and hypoallergenic.

- Flexipalette can be sterilized in autoclave at 135°C
- Flexipalette is absolutely opaque and reflection-free
- Flexipalette is gentle and soft in the patient's mouth
- Flexipalette can be folded in a 30° to 45° degree angle

Flexipalette A



Flexipalette B



Flexipalette C



Flexipalette Color Match



With Flexipalette

Without Flexipalette

A must have for
intra-oral photography!



Description	Item no.
Flexipalette, set of 4 forms	FP-3100
Flexipalette, form A	FP-3001
Flexipalette, form B	FP-3002
Flexipalette, form C	FP-3003
Flexipalette, Color Match	FP-3004




Discover an instructive video-clip demonstrating you how to use Flexipalette and its great advantages.

PHOTOCAD

A NEW GENERATION OF RETRACTORS



PC-3001
Photo-CAD LOW

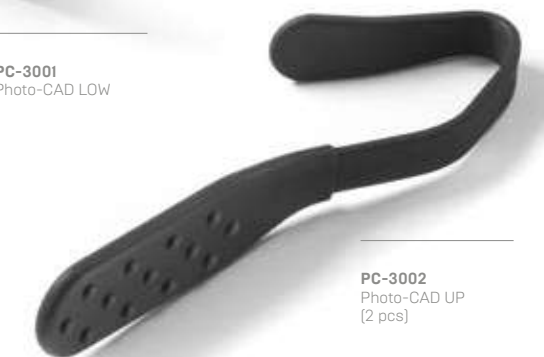

PC-3002
Photo-CAD UP
(2 pcs)

The innovative design of the Photo-CAD retractors allows for the best results with photography of the 4 quadrants and during digital impressions with Cad CAM systems. Photo-CAD can be ergonomically folded and adapted to the patient's mouth and to the needs of the photographer. They have a soft surface for the patients' best comfort. **The material is opaque, reflection-free** and can be sterilized.

Photo-CAD LOW is a totally unique tool which allows for the retraction of both the cheek and the tongue in the same use during an "occlusal" image capture of a lower quadrant. It allows to make the arch free by keeping aside the soft tissues. Photo-CAD LOW also makes optical scanned-impressions of lower preparations easier.

Photo-CAD UP is a cheek retractor which can be used every day. Its function when in use for dental photography will be to retract the vestibular soft tissues in the area upper pre-molar/ molar during "occlusal" image capture. It will be of great help too during optical scanned-impressions for Cad-CAM systems.

Description	Item no.
Photo-CAD, set of 3 pcs	PC-3100

SIDEKICK

IMPLANT DRIVERS CLINIC



**Nobel Biocare¹
Compatible driver**
D-400-NOB-S, 8 mm
D-400-NOB-M, 13 mm
D-400-NOB-L, 25 mm



**Straumann¹
Compatible driver**
D-400-STR-S, 8 mm
D-400-STR-M, 13 mm
D-400-STR-L, 25 mm



**Astra Tech¹
Compatible driver**
D-400-AST-S, 8 mm
D-400-AST-M, 13 mm
D-400-AST-L, 25 mm

Sidekick's box of screwdrivers for implantology gathers 9 models (each available in 3 lengths) that are compatible with the most popular implant systems on the market. The aim is that you have available, **anytime and immediately THE correct model and length of driver needed.**

Thanks to the wide choice of proposed geometries, in 99% of the cases you'll find a compatible model with other brands of implants' screws. A smart clipping mechanism allows you to quickly assemble/disassemble the driver's tip and the (very compact) head.

For completing the concept, Smile Line has created a special clamp forceps which tip is perfectly designed to match the profile of the driver's head for extra tightening or difficult un-tightening. Its use is extremely comfortable **with a direct access to the head's side** without removing the driver from its place.

Smile Line's box of implant drivers is available in two executions: the "basic" set that includes one driver tip of each in the medium size [9 tips + 1 head] or the "master" set that includes one driver tip of each in the three sizes [27 tips + 1 head].

Should you prefer to compose your own customized box, you may order all the components individually.



**Friadent¹
Compatible driver**
D-400-FRI-S, 8 mm
D-400-FRI-M, 13 mm
D-400-FRI-L, 25 mm



**3i¹
Compatible driver**
D-400-iii-S, 8 mm
D-400-iii-M, 13 mm
D-400-iii-L, 25 mm



**Camlog¹
Compatible driver**
D-400-CAM-S, 8 mm
D-400-CAM-M, 13 mm
D-400-CAM-L, 25 mm



**Ankylos¹
Compatible driver**
D-400-ANK-S, 8 mm
D-400-ANK-M, 13 mm
D-400-ANK-L, 25 mm



**SPI¹
Compatible driver**
D-400-SPI-S, 8 mm
D-400-SPI-M, 13 mm
D-400-SPI-L, 25 mm



Flat driver 1.4 mm
D-400-1.4-S, 8 mm
D-400-1.4-M, 13 mm
D-400-1.4-L, 25 mm



D-401
Driver's head



Ergonomic, universal and comfortable.



Description	Item no.
Sidekick implant drivers set, inclu. 9 tips 13 mm, 1 head, aluminium box	D-430-Basic
Sidekick implant drivers set, inclu. 27 tips [8-13-25 mm], 1 head, aluminium box	D-430-Master
Sidekick aluminium box, empty	D-499
Replacement o'ring for driver's head, 5 pcs	D-403

SIDEKICK CLAMP FORCEPS

The design of its jaws allows lateral access and grip of the screwdrivers ! In short, much better comfort and ergonomics!



D-402
Clamp forceps



DYNAMOMETRIC KEY

Adjustable: 15 / 20 / 25 / 35 Ncm

scale 1:1



D-490
Dynamometric key

Thanks to its adapter (#D-491), this dynamometric key is compatible with all Sidekick's drivers tips.

Note : owners of an original Straumann key will only need the adapter D-491 for being able to use all the Sidekick tips in any circumstances.



D-491
Adapter for dynamometric key

¹Nobel Biocare, Straumann, Astra Tech, 3i, Camlog, SPI, Ankylos, Friadent are all registered Trademarks.

VENEER ME

Thanks to veneerME, you'll be able to carry and ship and also manipulate your hyper-delicate veneers in the highest conditions of security!

Maximum security for your veneers.



Description	Item no.
veneerME, set of 3 complete boxes, 3 spare PU pads and 3 spare membrane-frames	2600
veneerME, spare PU pads (10 pcs)	2610
veneerME, spare membrane-frames (10 pcs)	2620
veneerME, Mini-set of 1 complete box, 1 spare PU gel pad and 1 spare membrane-frame	2600-1

SMILELINE 2023


THE STICKY POLYURETHANE GEL-PAD

The base of veneerME box is covered with a sticky pad out of polyurethane gel. Veneers are placed on this pad on their buccal surface, and adhere naturally. All the preparation process of the internal part of the veneer (acid etching, salinization) can be done in this position **without any fear of a veneer escaping or falling on the floor**. The pad's material is super resistant and can be re-used several times (it is easily replaceable in case of deterioration).

THE SECURITY MEMBRANE-FRAME

Once the veneers are secured on the pad, the membrane-frame is clipped on the base of veneerME box for extra fixation. The xx-thin membrane, made out of nylon, delicately adapts to the shape of the veneers. This done, you just need to close the lid and your master work is ready for delivery. The membrane-frame can be replaced -if needed- after having been used several times.

HUGE ADVANTAGE: THE MEMBRANE IS PERMEABLE! In other words the clinician can proceed with rinsing and drying after having replaced the membrane-frame on the base. This way the clinician can work in the best conditions and in total security.

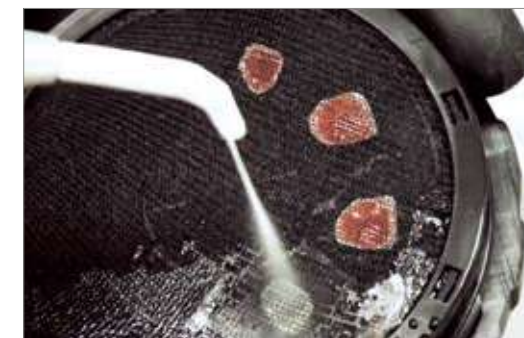

HOLDING GEL

Gel for temporarily fixation.

A small drop on the die is enough for holding very well – and temporarily – your delicate porcelain veneers as well as your full-ceramic restorations on the model.

This gel is essential for **checking the contact points, for the application of correction powders with great ease and for restorations that oppose each other.**

Transparent, so it doesn't alter the shade – it completely dissolves in water and doesn't leave any residue after firing.



Description	Item no.
Holding gel, 2 x 6 ml	5110



MY SHADE GUIDE FOR COMPOSITES

EASY MAKING OF DENTIN ENAMEL SAMPLES

Thanks to the special design of the two-parts mould, you have the possibility to produce very easily and in a professional way your **composite teeth samples dentin-enamel**.

A small portion of enamel is placed in the first mould, then pressed with help of the second mould which is in fact presenting a positive shape of dentin. Result after light-curing of this first step will be a composite veneer out of enamel material. The second step will need the filling of the hollow part of the veneer with the selected shade of dentinal composite material. Simple, fast and professional!

Two sets are proposed (Minikit and Masterkit), offering both the needed accessories for the production of your customized shade guide.

SMILELINE 2023



SAMPLER

This very precise tool allows you to produce composite tabs from your personal materials in a very simple and professional way. Thanks to the graduation, you can decide of the thickness of your samples.


Should you wish, black tablets out of laminated card stock exist and give you the possibility of gluing and indexing all your tabs for complete customization.

Sampler for composite tabs.



Transparencies according to different thicknesses

Description	N°Art.
My Shade Guide, Mini-kit: 1 holder, 5 sticks&clips, 1 double mould in silicone	7600-N
My Shade Guide, Master-kit: solid aluminum base with PC lid, 5 holders, 25 sticks&clips, 1 double mould in silicone	7500-Master
Double mould in silicone	7630
Sticks & clips, 100 pcs	7520-N
Sticks & clips (25 pcs) + holders (5 pcs)	7522-N
Clips, 100 pcs	7525-N

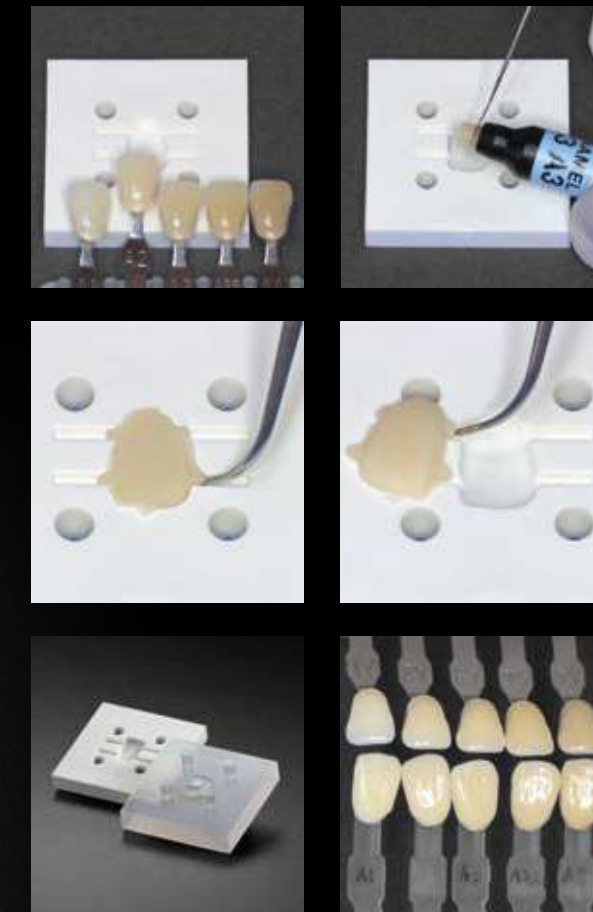

Description	Item no.
Sampler only	7015
Indexing tablet, 5 pcs	7016

COMPO-VIBES

SMART VIBRATIONS !



The guarantee of an easier, faster and better composite's modeling ! Thanks to the micro-vibrations it generates, Compo-Vibes optimises the material's thixotropic properties and makes your work easier.

- Better adaptation of the composite to the cavity wall, no detachment, less air bubbles
- Composite much less sticky than with a traditional instrument
- Modeling and sculpture of composites much easier and fast
- Increased control of material thicknesses
- More precise application

Description	Item no.
Compo-Vibes Basic Set: Compo-Vibes unit with one instrument "Tip medium"	03000-CV-Basic
Compo-Vibes Master Set: Compo-Vibes unit with 4 instruments: "Spatula, Tip large, Tip medium, Tip thin"	03000-CV-Master
Set of 4 instruments : "Spatula, Tip large, Tip medium, Tip thin"	03010-CV
Spatula, 2 pcs	03011-CV
Instrument "Tip large", 2 pcs	03012-CV
Instrument "Tip medium", 2 pcs	03013-CV
Instrument "Tip thin", 2 pcs	03014-CV
Compo Brush, spare flat tips for Compo-Vibes #Anterieur (30pcs)	4651-CV-ANT
Battery AAAA 1,5 V	01030-MV
Disposable protective sleeves, 200 pcs	3011-KL

SMILELINE 2023



COMPO-BRUSH

COMPOSITE STAINING

The Compo-Brushes developed with Styleitaliano are intended for composite staining.

Their black synthetic fiber allow for **extreme precision with the best contrast**. A working quality superior to standard brushes you find normally in the composite kits. Brush tips hold by friction on their handle.

#000 Posterior brush tip can easily be bent by the user for a better access to the posterior areas.

Compo-Brush set contents two complete brushes and two replacement brush-tips of each model. Single brushes as well as brush tips refills are available separately.

Description	Item no.
Compo-Brush, set of 2 brushes and 2 spare tips ea. (2 x ANT & 2 x POST)	4650-Set
Compo-Brush, tips #000 for posteriors (50 pcs)	4651-POST
Compo-Brush, flat tips for anteriors (50 pcs)	4651-ANT




4650-POST
Compo-Brush
#000, complete

4650-ANT
Compo-Brush
flat, complete



ORGAN.EZ

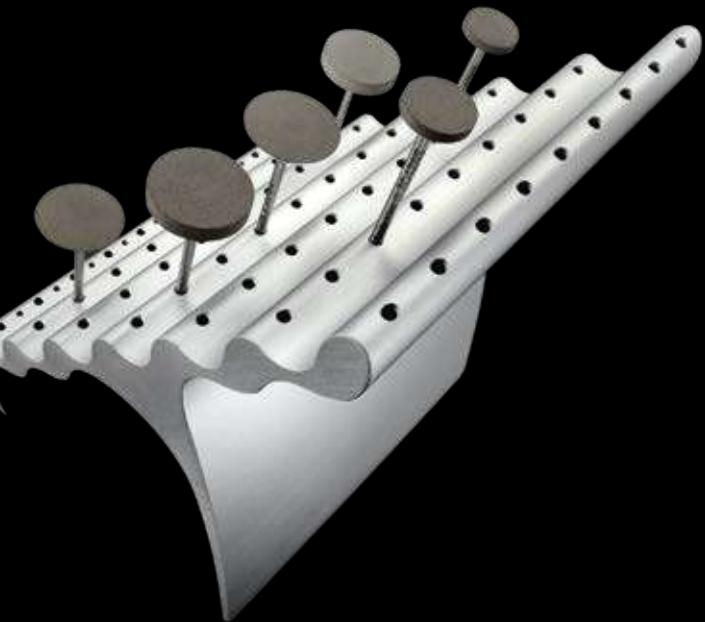

Solid aluminium burstand.



At first sight, some will see a sculpture placed on the bench. Don't be misled: ORGAN.EZ is first and foremost extremely useful and ergonomic.

ORGAN.EZ allows you to **systematically arrange** your burs, diamonds, polishers with a simple notion – knowing instantly where to find them instead of having to search.

Made out of **solid aluminium**, ORGAN.EZ is very stable where you place it. For a maximum comfort when taking or replacing the burs, the entrance of each hole starts with a slight angle.

A maximum of 63 instruments with a shaft 2.35 can be placed. On the first level, you'll notice 8 additional holes intended for the placement of your air-turbine burs.

8100
Organ.ez, solid aluminium stand for HP & FG burs

SMILELINE 2019

SILICONE MOULDS

N-ARCHITECT ANTERIORS



Made out of transparent silicone, these moulds allow you to easily and rapidly produce your own wax or composite veneers.

The **ANTERIOR** assortment is a collection of beautiful basic shapes: triangular, oval, rectangular and includes a universal lower anterior form, all available in two sizes for a total of 8 moulds in the set.





2111
New Architect, set of 8 anterior silicone moulds

GOLDEN SECTION DIVIDER

The golden proportion is a ratio which result is called mathematically Phi (1.618) or phi (0.618). It is proven that **the universe uses this proportion for its own construction**, from the genetic sequences up to the level of galaxies, passing by construction of flowers, insects, a hand, a face or the whole body.


As we are able to divide a line by its middle or draw a perpendicular, we all recognize the Golden Proportion as a constant for harmony and balance.

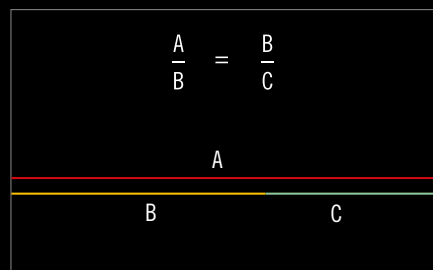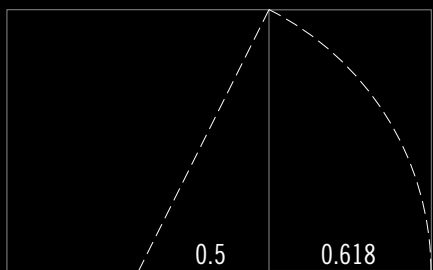

Description	Item no.
Golden Section Divider set, including an assortment of 6 grids blocs, instructions	6010

For the determination of the Golden Proportion.

1 1 2 3 5 8 13 21 34 55 89 144 233 377 610 987 1597 2584 4181 6765 10946 17711 28657 46368 75025 121393 196418 317811 514204 832017 1346221 2178232 3524453 5702684 9229837 14932511 24155348 39088165 63243486 102381835 165625221 268007106 433688941 701714166 1135402067 1837116233 2972518300 4809614533 7781632833 12590170146 20371787679 32961957825 53333735404 86305523083 139669260487 225984783570 365651443657 591636227227 957616010800 1549253234457 2506869245257 4056116479657 6562985714914 10619096954571 17185213674185 27804310629099 44989523583674 72793834202773 117603146881872 190396974028645 308000120910417 498396294939292 806396415850709 1304792610789901 2113188805729193 3418580416518802 5531773027208000 8954961832927193 14486735949145195 23441708006353198 37928674155498393 61370412162641591 100000000000000000

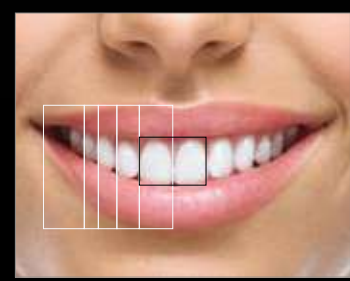

FASCINATING TOOL
Golden Section Divider will **allow you to easily assimilate "physically" the concept**. Slightly or widely open, the Golden Section Divider indicates always two segments respecting the golden proportion. Made out of the highest quality stainless steel and very carefully assembled by hand, the divider is a wonder of precision mechanics. It is supplied together with an assortment of the very useful blocs of grids and of course with instructions of use.



THE GOLDEN PROPORTION IN DENTISTRY
Many studies mention the golden proportion as being one of the ingredients for the harmony of a smile and dentition (see images). The Golden Section Divider and the grids allow you to discover in a few seconds the theoretical ideal goal in comparison to the real situation in the mouth. These information are very exciting and very useful for:

- aesthetic anterior cases
- dentures
- positioning-angle-axis of teeth
- gingivectomy
- wax-up control



The "Golden Section Grids" printed on paper (assortment of blocs supplied in the kit) are very useful for registrations in the mouth of a patient as well as when working with a model.

STYLEITALIANO

MORE THAN JUST AN INSPIRATION FOR US

The study group Styleitaliano started with a meeting between two friends, Dr. Walter Devoto and Prof. Angelo Putignano, both worldwide renowned experts in restorative aesthetic dentistry. Their goal is to share with the most people simple and proven techniques allowing everyone to improve and tend towards success.

Since then, the study group met a huge development for becoming one of the most respected and followed group internationally.

WWW.STYLEITALIANO.ORG

WWW.FACEBOOK.COM/STYLEITALIANO

Smile Line and Styleitaliano met in 2011. While we were focusing on the development of tools and accessories dedicated to the dental laboratory, Styleitaliano allowed us -thanks to their innovative and brilliant ideas- to open the door of dental offices.

We wish to express here our warm thanks to every person in the Styleitaliano group for their kindness, availability and humility.





Smile Line ☺

made in switzerland

Smile Line SA
Z.I. La Clef 5
PO Box 229
CH-2610 St-Imier · Switzerland
Tel. +41 32 944 22 44
Fax +41 32 944 22 45
mail@smileline.ch
www.smileline.ch