# Porcelain Restorations

Dental - Labor | Dental Laboratory | Laboratoire Dentaire | Laboratorio Dentale | Laboratorio Dental  
Zahnarzt - Praxis | Dental Practice | Cabinet Dentaire | Studio Dentistico | Clínica Dental

edenta

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All rights reserved.

Reproduction, also by extract and reproducing photos are only permitted with written authorization of EDENTA AG, Switzerland. All offers, orders and deliveries are subject to Edenta's „General sales and Delivery Terms“. We reserve the right to modify our range of products and their design as well as to deviate from the illustrations and data shown.

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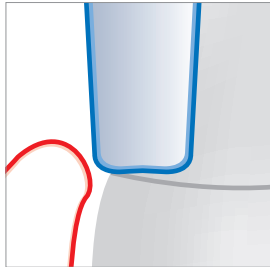
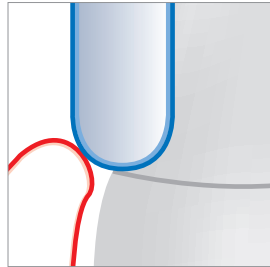
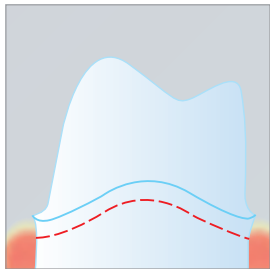
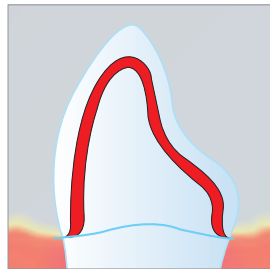
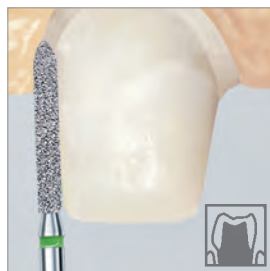
Rotary instruments as integrated components in the dental practice for maximum safety during preparation and optimum aesthetics. Also specially designed for all-porcelain restorations with zircon oxide.

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**Preparing teeth for crowns**

The tooth is generally prepared according to the usual principles for shoulder and deep chamfer preparations with a preparation depth of 1 mm at the preparation margin.

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Specially developed and patented WhiteTiger diamonds, enable a controlled and minimally invasive reduction of the tooth structure. For efficient preparation of porcelain restorations with either a shoulder or deep chamfer preparation in the anterior and posterior regions.

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## WhiteTIGER Diamonds for crown preparation

The multilayered diamond coating and patented LLD bonder ensure increased material reduction with excellent chip flow. There is greater reduction of tooth structure compared with standard diamonds, saving time during treatment as well as reducing patient discomfort.

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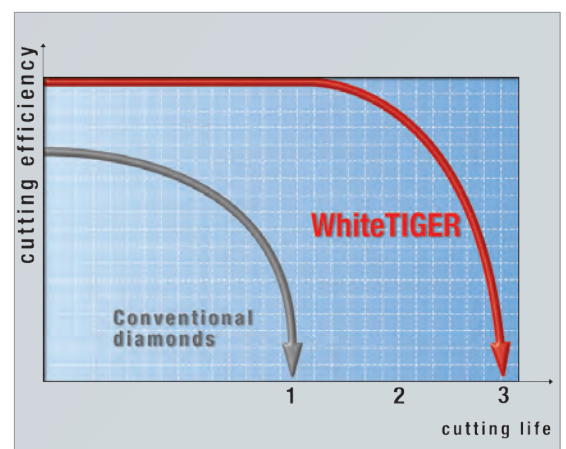
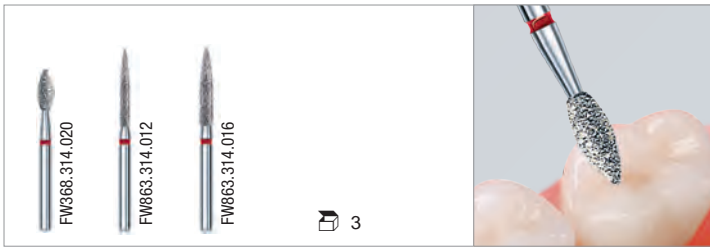
### Occlusal Adjustment and Polishing

Fine-grit White Tiger diamond instruments are used for adjusting occlusal contact points. Ceramic polishers are used for polishing.

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### Trepanation with WhiteTiger diamonds

An access cavity is prepared in two stages. ① The facing porcelain is first removed with a diamond instrument without perforating the porcelain framework. ② The porcelain framework is then perforated with a diamond instrument leaving a gap of 0.5 mm from the bonded porcelain (to prevent bonded porcelain breaking off).

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### Polishing system for porcelain restorations

3 stage polishing system for top quality pre-polishing, high-lustre and super high-lustre polishing of all-porcelain restorations. The specially developed polisher matrix guarantees an optimum surface sheen without the use of polishing paste thereby ruling out the need for an additional glaze firing after adjusting the porcelain.

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3 - step diamond polishing system suitable for all porcelain fillings, partial and full-crown restorations.

- Step 1 = Green:** Finishing, for rapid material reduction
- Step 2 = Blue:** Polishing, for smoothing the surface
- Step 3 = Yellow:** High - lustre polishing

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3 - step diamond polishing system suitable for all porcelain fillings, partial and full-crown restorations. Polishers with natural diamond grit suitable for finishing, polishing and high-lustre polishing of porcelain surfaces. A second glaze firing is not necessary following adjustment.

- Step 1 = Green:** Finishing, for rapid material reduction.
- Step 2 = Blue:** Polishing, for smoothing the surface
- Step 3 = Yellow:** High-lustre polishing

**Application: Use coolant spray when polishing.**

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2 - stage finishing and polishing system for all types of porcelain fillings, partial and full crown restorations. Polishers for finishing and polishing porcelain surfaces. Further glaze firing is unnecessary after adjustment.

Porcelain polisher **light grey** = abrasive prepolish, smoothes surfaces ground with diamonds.

Porcelain polisher **pink** = slightly abrasive polish, maintains the anatomical structure.

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Polishers for all porcelain fillings, partial and full-crown restorations. For polishing porcelain surfaces to a high-lustre, producing a natural sheen without the use of polishing paste. A second glaze firing is not necessary following adjustment.

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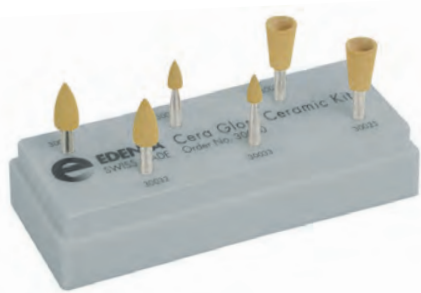
Set enthält 6 CeraGloss RA Polierer.

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Unique diamond - impregnated surface treatment system for high - gloss polishing of surgery ceramics. Time - saving, convenient single step process.

Natural, high temperature - resistant rubber with top - grade natural diamond particles which ensure rapid surface reduction combined with a high - gloss result of unsurpassed excellence.

Set contains 6 CeraGloss RA polishers.

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Kit consists of 9 assorted instruments in RA and FG shank, specifically designed for the Porcelain Veneer Technique. The fine diamond (red ring) is used in preparation of the tooth and the special silicone CeraPol Plus diamond impregnated points are used to obtain a smooth and satin finish.

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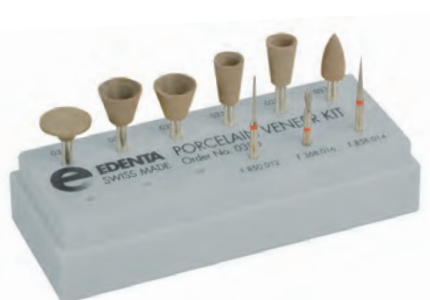
### Removing an all-porcelain restoration (crown cutter)

To remove an all-porcelain restoration, cut the restoration open along the length of the axial wall to the centre of the occlusal surface or incisal edge using the sinter diamond disc order no.5114RA. Then bend the restoration open with a suitable instrument to make it fracture. No. 9920 Cleaning stone

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### Removing a low-fusing porcelain restoration.

Crown cutter C4AK for all types of metal alloys, titanium and low-fusing porcelain and composite facings: **C4AK Crown cutter with a new MultiCut design.** Separates crowns and bridges easily and quickly in one working stage. These crown cutters are not suited for cutting all - ceramics restorations of extremely hard ceramics, such as zirconium - oxide ceramics.

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### Finishing an unsintered bridge framework

Smoothing the connectors after milling the bridge framework in a CAD milling unit with tungsten carbide cutters and polishers CeraGloss

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### Optimum material reduction and excellent surface quality

when preparing porcelain and composite facings. Minimum vibration when preparing results in optimum surface quality, superior to preparing with diamond instruments. Controlled material reduction and increased cost-effectiveness.

- opt. 25.000 **Verblendkunststoffe**, Facing composites, Matériaux cosmétiques  
Resine per faccette composito, Resinas de recubrimiento
- opt. 15.000 **Niedrig schmelzende Keramik vor dem Glanzbrand**.  
Low-fusing porcelain before glazing.  
Céramique à basse fusion avant cuisson de glaçage.  
Ceramica a basso punto di fusione prima della glasatura.  
Cerámica de baja fusión antes de la cocción de glaseado.

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
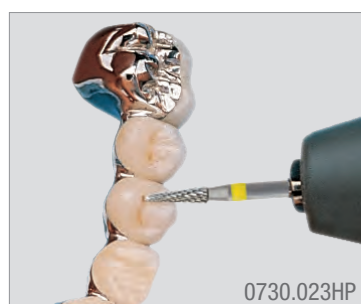
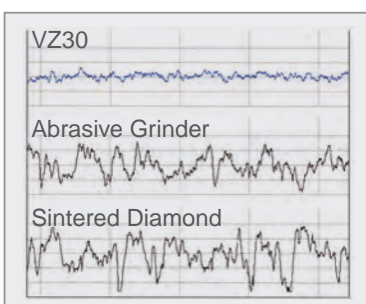
Diamond burs for trimming sintered zircon oxide

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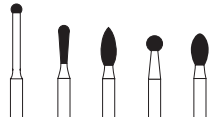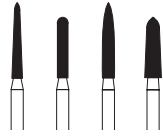
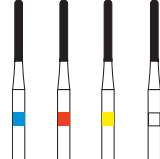
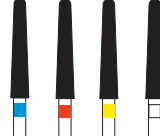
## Fitting the sintered zircon oxide framework

By trimming the high spots with diamond instruments.

For trimming and adjusting sintered /hot isostatically pressed (HIP) zircon oxide using a laboratory turbine and cooling with water.

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## MT Milling Diamonds

For use in a milling machine with air turbine and water cooling when milling all-ceramic primary telescopes, telescopic implant components, all-ceramic abutments and zircon oxide primary copings with optimum surface quality.

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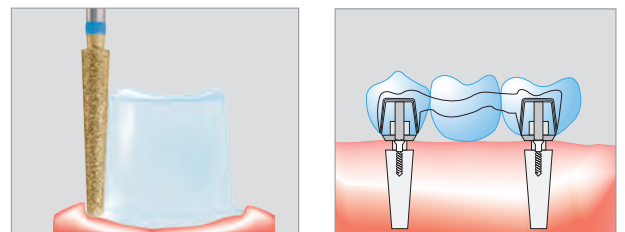
Instruments for ZrO<sub>2</sub>, porcelain, all-porcelain and pressable porcelain restorations

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#### Rough trimming

Sintered diamond with organic bonder, for grinding of ceramics and metals, enables rough trimming without heat build-up in the restoration. Suitable for rough trimming and pressable porcelain sprues.

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#### Initial contouring

Abrasive trimmer with natural diamond grit. Ideal for initial contouring and preparing without heat-build in the restoration. Suitable for preparing larger vestibular and oral surfaces.

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### Polish

3-step diamond-polishing system for working on all types of ceramic without heat generation. A new binding component as well as the natural diamonds allow for optimum polishing results, achieved without any polishing paste.

1. step: Green = coarse grit
2. step: Blue = medium-coarse grit
3. step: Yellow = super - fine grit

for pre - grinding, for quick and bulk material reduction.  
for smoothing the surfaces and preparation for final polishing.  
for high - shine polishing without additional final glaze.

### Instructions for use:

- Apply only light pressure when polishing.
- In order to achieve optimal polishing results, it is important to use the three polishers in the sequence as mentioned.

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
**Step 1 = Blue – Coarse grit**  
for pre-grinding, for quick and bulk material reduction.

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
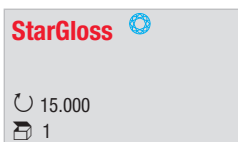
**Step 2 = Pink – Medium grit**  
for smoothing the surfaces and preparation for final glaze.

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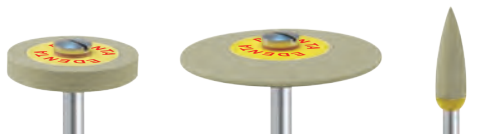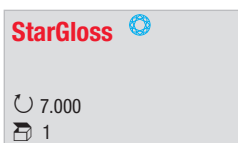**Step 3 = Grey – Superfine grit**  
for high-shine polishing, without using polishing paste or needing an extra glaze firing.

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### Polish

New ceramic systems and their processing techniques require new polishing systems. Specifically developed for gently polishing porcelain and all-ceramic restorations (zirconia).

Natural diamond grit with a highly condensed, long-mesh synthetic bonder, produce an outstanding surface texture. Finish, smooth and high-lustre polish without using polishing paste or needing an extra glaze firing.

**These StarGloss polishers are guaranteed to polish all-ceramics gently and successfully.**

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2- step finishing and polishing system suitable for all partial porcelain and full-crown restorations.

**Step 1 = grey - white:** abrasive, eliminates scratches and smoothes the surface.

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2- step finishing and polishing system suitable for all partial porcelain and full-crown restorations.

**Step 2 = pink:** retains the structure and provides a final shine.

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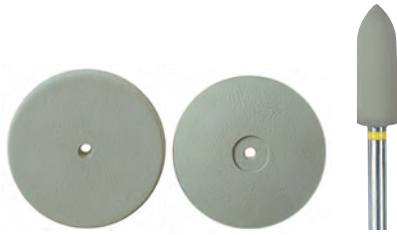
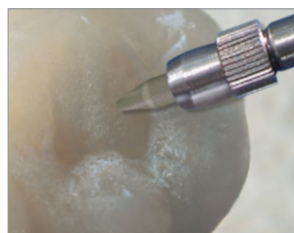
Polishers for all partial porcelain and full-crown restorations. For polishing porcelain surfaces to a high-lustre, producing a natural sheen without the use of polishing paste.

A second glaze firing is not necessary following adjustment.

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Polisher with diamond impregnation for preparation of occlusal surfaces on ceramics.

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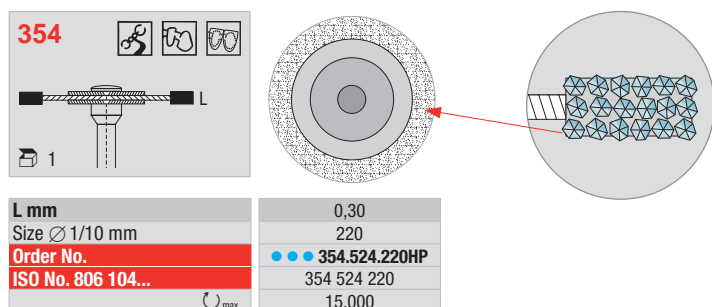
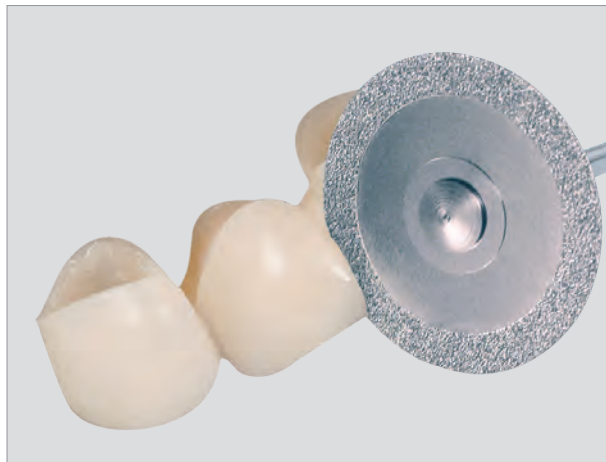
## Properties

- Multilayered diamond coating for grinding edges and large areas
- High removal and cutting capacity • Prepare without striations
- Flexible, highly stable disc • Long service life

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Diamond disc with a multilayered, electroformed diamond coating on the rim. Manufactured by electroforming.

The diamond coating on the rim ensures a cutting performance similar to that of a sintered diamond disc, but with a higher degree of flexibility.

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## Applications

Used in the dental laboratory for Cutting off, Pretrimming, Separating and Contouring of Porcelain, Zirconia, Precious metal as well as CrCo alloys and Composites.

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**Zum Abtragen 20.000 upm / Zum Polieren 10.000 upm**

Set consists of 1 dozen grey-white wheels for pre-polishing.

1 dozen pink wheels for final polishing and 1 dozen grey wheels for high-gloss polishing of porcelain, 2 TC cutters, cross cut superfine and 3 mandrels.

Important note: rpm recommendations

**i.e. for grinding and abrasive work 20.000 rpm / for final polish – below 10.000 rpm**

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Unique diamond-impregnated surface treatment system for contouring & high-gloss polishing of laboratory ceramics.

Time-saving, convenient two-step process.

Synthetic, high temperature-resistant rubber with a unique concentration of top-grade natural diamond particles which ensure rapid surface reduction combined with a high-gloss result of unsurpassed excellence.

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Special rotary instruments for preparing porcelain, pressable porcelain and all-porcelain without heat build-up to prevent microcracks. Preparation is completed in three stages:

- |                              |                             |
|------------------------------|-----------------------------|
| <b>1 Initial preparation</b> | SuperMax rotary instruments |
| <b>2 Precontouring</b>       | CeraPro rotary instruments  |
| <b>3 Polishing</b>           | CeraGloss polishers         |

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**Oberflächenschliff:** HM-Fräser

**Politur:** CeraGloss Diamantpolierer

22 instruments for preparing porcelain, all-porcelain and pressable porcelain restorations.

**Initial preparation:** SuperMax,

**Precontouring:** Diamond rotary instruments

**Surface preparation:** Tungsten carbide cutters

**Polishing:** CeraGloss diamond polisher

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